

High Speed Simultaneous Analysis of Amino Acids in Foods Using Automatic Pretreatment Function

T. Yoshioka

User Benefits

- ◆ The analysis time is significantly reduced compared to the post-column derivatization method.
- ◆ The analytical method is very simple and easy since the burdensome derivatization process can be performed automatically.
- ◆ Users can analyze amino acids with a simple HPLC system.

Introduction

Although the post-column derivatization method has been commonly used for amino acid analysis with HPLC, its disadvantages include the long analysis time due to the characteristics of the column and high cost from the complex instrument configuration. On the other hand, although pre-column derivatization enables fast analysis with simple instrument configuration, its problems include the burdensome derivatization operation and the effect of the sample matrix.

In Application News 01-00007, we introduced a pre-column derivatization method for amino acid analysis using the automatic pretreatment function. With this feature, the burdensome derivatization process can be performed automatically. In this article, we introduce the examples of analyzing amino acids in various foods using the method with automatic function.

Automatic Pre-Column Derivatization

Nexera XR is equipped with an automatic pretreatment function that enables users to configure desired operations including sample dilution and reagent addition. For this study, we set the system to automatically mix the sample and derivatization reagent in the autosampler needle. Fig. 1 shows the flow of the derivatization, and Table 1 the preparation method for the derivatization reagents. For detailed pretreatment program parameters, please refer to Application News 01-00007.

Analysis of Mixed Standard Amino Acid Solution

Fig. 2 shows the chromatogram of a mixed standard solution of 20 proteinogenic amino acids. The 20 components could be separated in about 14 minutes. Tables 2 to 4 on the next page show the analytical conditions.

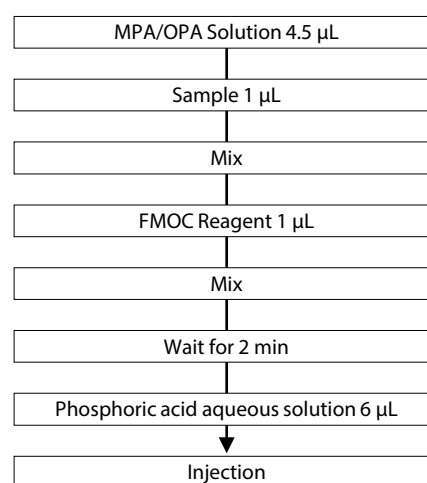


Fig. 1 Derivatization Flow with the Automatic Pretreatment Function

Table 1 Preparation of Derivatization Reagents

- 0.1 mol/L Borate buffer
Add 0.62 g of boric acid and 0.20 g of sodium hydroxide into 100 mL of ultrapure water.
- Mercaptopropionic acid Reagent (MPA Reagent)
Add 10 µL of 3-mercaptopropionic acid into 10 mL of 0.1 mol/L borate buffer.
- OPA Reagent
Add 0.3 mL of ethanol into 10 mg of o-phthalaldehyde and dissolve completely. Then add 0.7 mL of 0.1 mol/L borate buffer and 4 mL of ultrapure water.
- MPA / OPA Solution
Mix 600 µL of MPA Reagent and 300 µL OPA Reagent.
- FMOc Reagent
Dissolve 10 mg of 9-fluorenylmethyl chloroformate into 100 mL of acetonitrile.
- Phosphoric acid aqueous solution
Add 0.5 mL of phosphoric acid into 100 mL of pure water.

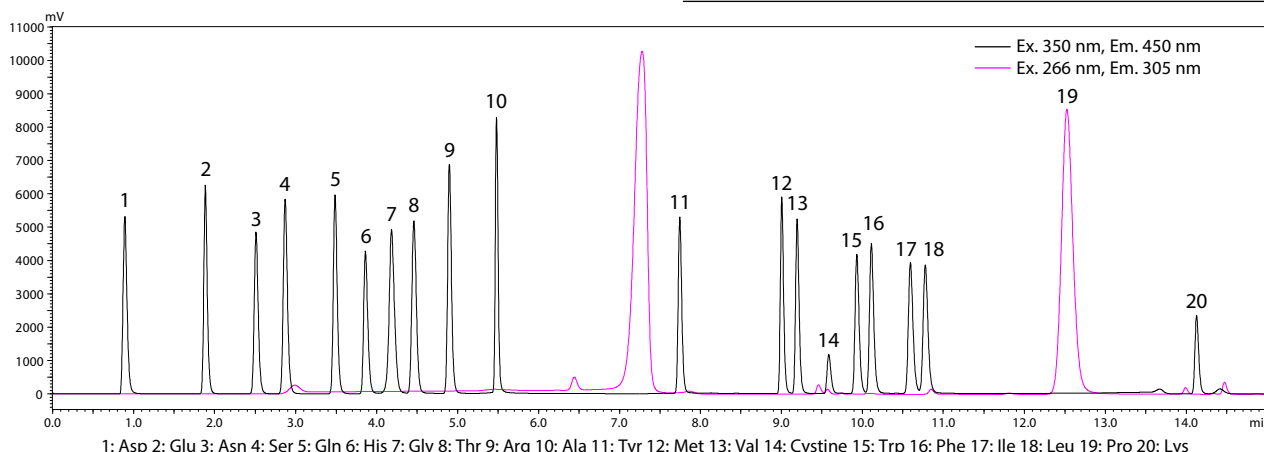


Fig. 2 Simultaneous Analysis of 20 Proteinogenic Amino Acids (100 µmol/L each)

Table 2 Analytical Conditions

System	: Nexera XR
Column	: Shim-pack XR-ODS II ¹ 100 mm × 3.0 mm I.D., 2.2 μm
Mode	: Low pressure gradient
Mobile phase	: A) 20 mmol/L (Sodium) acetate buffer (pH 6) B) Water/Acetonitrile = 1 : 9 C) 20 mmol/L (Sodium) acetate buffer (pH 5) containing 0.5 mmol/L EDTA-2Na
Flow rate	: 1.0 mL/min
Column temp.	: 35 °C
Injection volume	: 1 μL
Sample cooler	: 4 °C
Detection	: Fluorescence detector (Cell temp. : 35 °C) Ch1) Ex. 350 nm, Em. 450 nm Ch2) Ex. 266 nm, Em. 305 nm

*1 P/N 228-41624-92

Table 3 Preparation of Mobile Phases

- Mobile Phase A
Add 2.67 g of sodium acetate trihydrate and 41 μL of acetic acid into 1000 mL of pure water.
- Mobile Phase C
Add 0.19 g of EDTA-2Na, 2.03 g of sodium acetate trihydrate and 308 μL of acetic acid into 1000 mL of pure water.

Table 4 Time Program

Time (min)	A.conc.	B.conc.	C.conc.
0	95	5	0
0.2	93	7	0
1	93	7	0
4	87	13	0
5	0	15	85
7.5	0	30	70
12	0	35	65
14	0	45	55
14.01	0	95	5
16.99	0	95	5
17	95	5	0
19.5	95	5	0

■ Analysis of Amino Acids in Foods

Amino acids are one of the key nutrients in food that contribute to a variety of physiological functions and also to the taste of food. We here introduce examples of the analysis of various foods (Figs. 3 to 12), and Figs. 13 to 16 show the pretreatment procedures for those analyses.

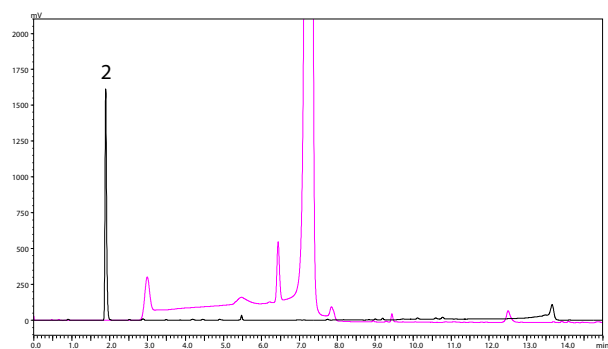


Fig. 3 Chromatogram of Consommé

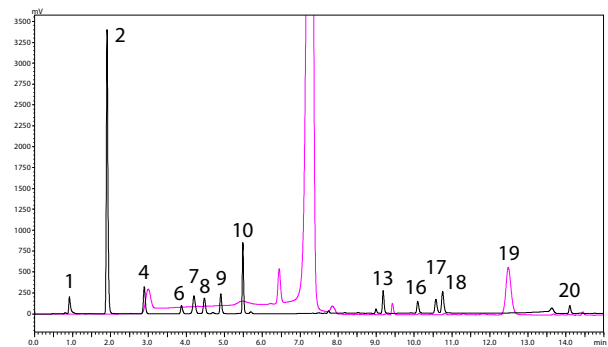


Fig. 4 Chromatogram of Japanese Soup Base

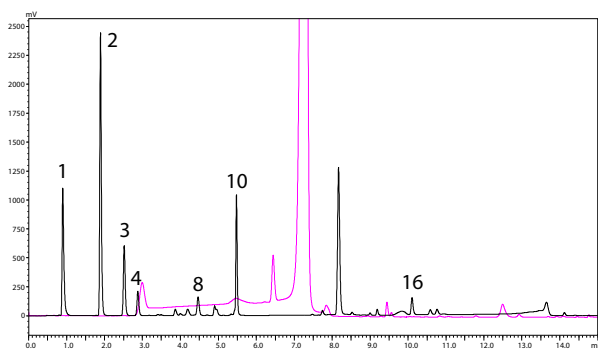


Fig. 5 Chromatogram of Ketchup

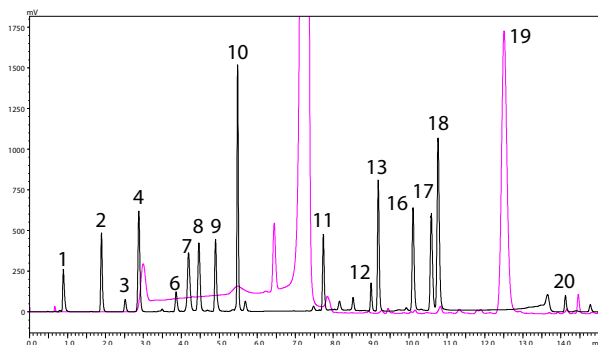


Fig. 6 Chromatogram of Miso

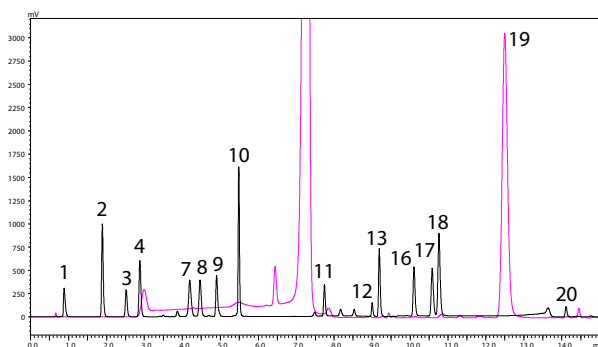


Fig. 7 Chromatogram of Gochujang

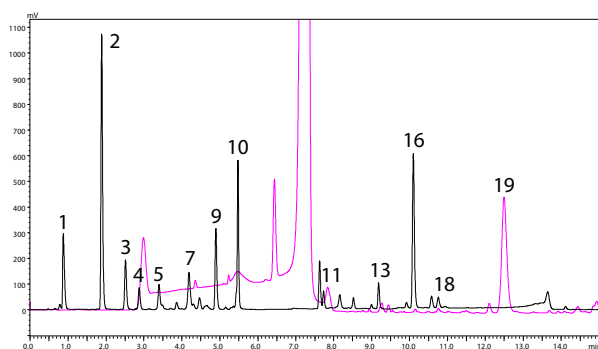


Fig. 8 Chromatogram of Peanut Butter

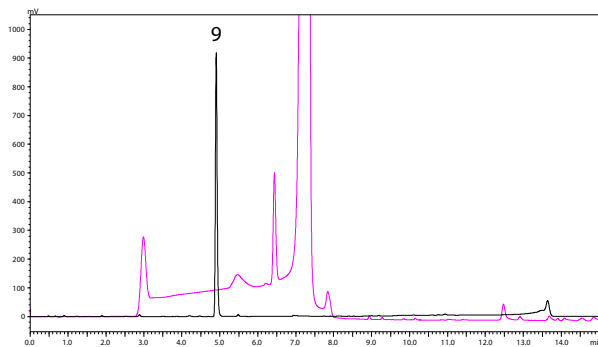


Fig. 9 Chromatogram of Energy Drink

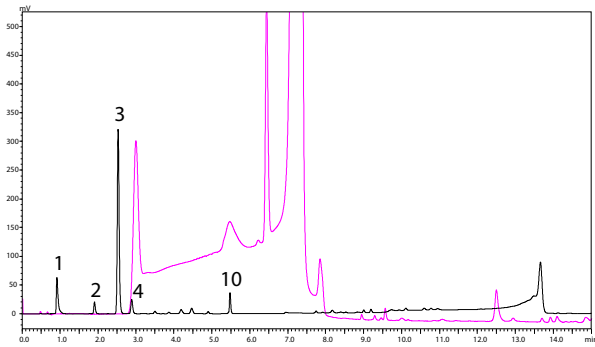


Fig. 10 Chromatogram of Peach Juice

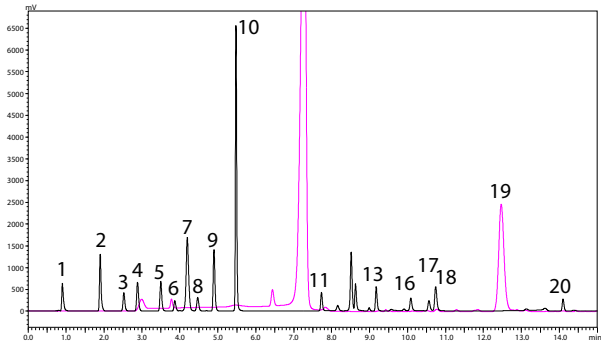


Fig. 11 Chromatogram of Sake

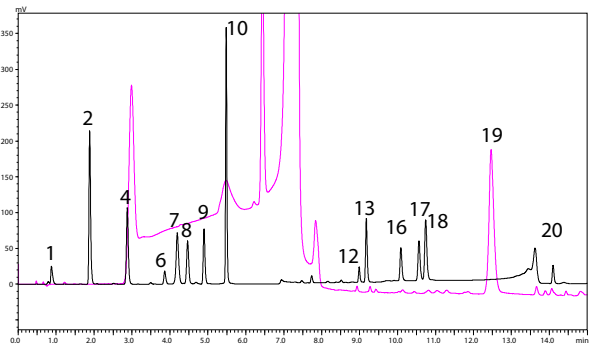


Fig. 12 Chromatogram of Soy Sauce

Sample 100 μ L

— Dilute with 1 mmol/L HCl up to 10 mL

Centrifugation at 3500 rpm for 10 min

Fig. 13 Pretreatment of Energy Drink, Peach Juice, and Sake

Sample 100 μ L

— Dilute with 1 mmol/L HCl up to 100 mL

Centrifugation at 3500 rpm for 10 min

Fig. 14 Pretreatment of Japanese Soup Base and Soy Sauce

Sample 1 g

— Dilute with water up to 10 mL

Vortex for 5 min

Centrifugation at 3500 rpm for 10 min

Supernatant 200 μ L

— Methanol 800 μ L

Vortex for 5 min

Centrifugation at 8000 rpm for 10 min

Supernatant 500 μ L

Evaporation to dryness

— 1 mmol/L HCl 500 μ L

Fig. 15 Pretreatment of Peanut Butter

Sample 1 g

— Dilute with water up to 10 mL

Vortex for 5 min

Centrifugation at 3500 rpm for 10 min

Supernatant 200 μ L

— Methanol 800 μ L

Vortex for 5 min

Centrifugation at 8000 rpm for 10 min

Supernatant 100 μ L

Evaporation to dryness

— 1 mmol/L HCl 500 μ L

Fig. 16 Pretreatment of Consommé, Ketchup, Miso, and Gochujang

■ Spike-and-Recovery Test

In general, it is assumed that the derivatization efficiency is affected by the sample matrix in the pre-column derivatization method. To check the effect of the sample matrix, we conducted spike-and-recovery tests on various samples to evaluate their recovery rates. Fig. 17 shows the results. Good recovery rates were obtained for many samples, indicating that the pre-column derivatization method can be applied to various foods.

■ Conclusion

The pre-column derivatization method introduced in Application News 01-00007 is a very simple analytical method that enables burdensome pretreatment process to be performed automatically. In this article, we introduced the examples of analyzing various foods using the analytical method. With this method, the amino acid analysis of foods can be conducted quickly and easily.

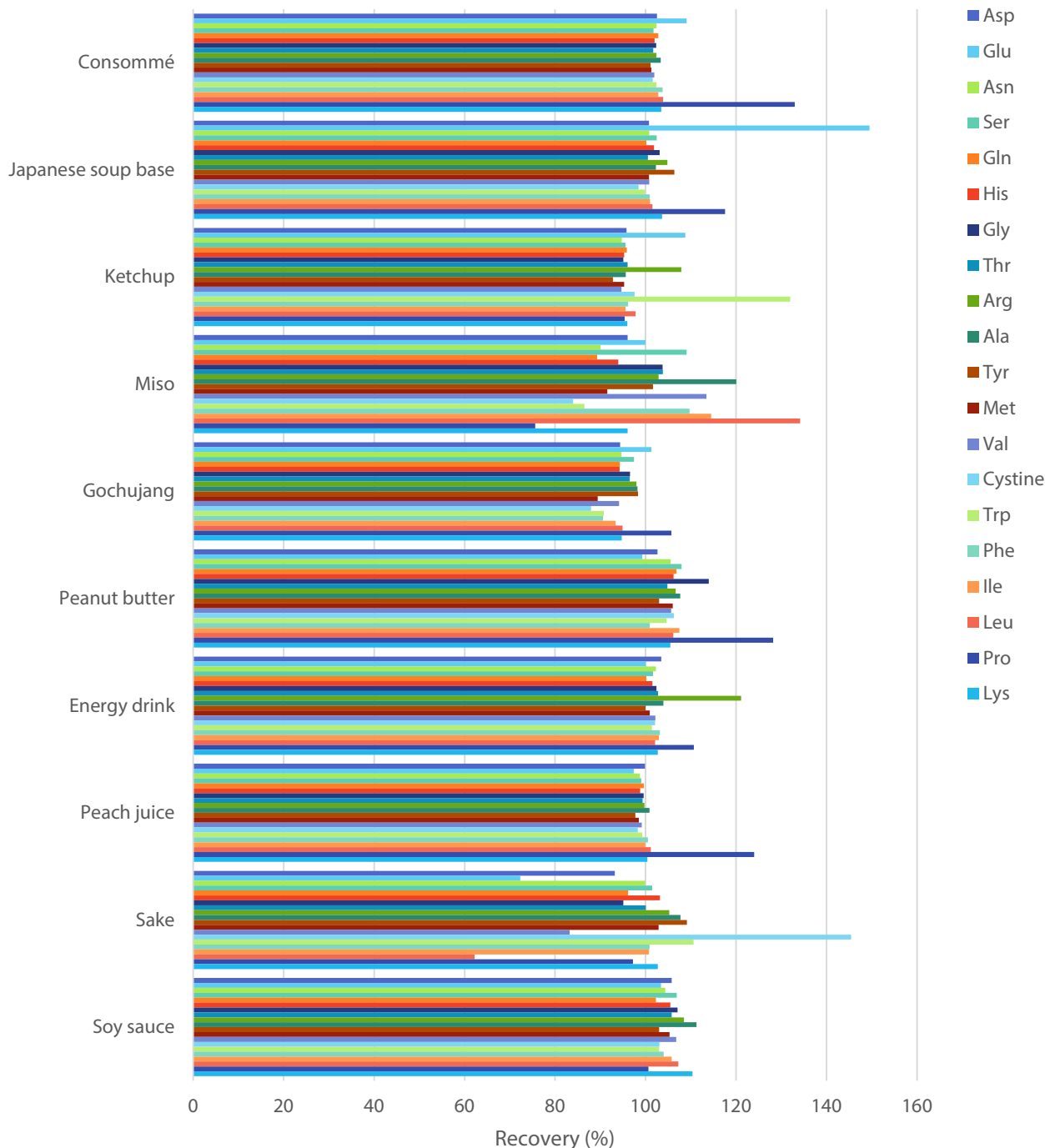


Fig. 17 Results of Spike-and-Recovery Tests

Nexera is a trademark of Shimadzu Corporation or its affiliated companies in Japan and/or other countries.



Shimadzu Corporation
Analytical & Measuring Instruments Division
Global Application Development Center

www.shimadzu.com/an/

For Research Use Only. Not for use in diagnostic procedure.

This publication may contain references to products that are not available in your country. Please contact us to check the availability of these products in your country.

The content of this publication shall not be reproduced, altered or sold for any commercial purpose without the written approval of Shimadzu. See <http://www.shimadzu.com/about/trademarks/index.html> for details.

Third party trademarks and trade names may be used in this publication to refer to either the entities or their products/services, whether or not they are used with trademark symbol "TM" or "®".

The information contained herein is provided to you "as is" without warranty of any kind including without limitation warranties as to its accuracy or completeness. Shimadzu does not assume any responsibility or liability for any damage, whether direct or indirect, relating to the use of this publication. This publication is based upon the information available to Shimadzu on or before the date of publication, and subject to change without notice.

01-00130-EN

First Edition: Sep. 2021